

Silicone Foam Dressing

Product information

Silicone foam dressing is an absorbent, atraumatic dressing made from absorbent core and silicone gel. The wound contact surface of Silicone foam dressing is coated with a layer of soft silicone gel that does not stick to the surface of a wound or cause trauma to delicate new tissue upon removal. This soft silicone layer is also slightly tacky, which facilitates application and retention of the dressing to intact skin, but does not cause epidermal stripping or pain on removal. The outer surface of the product is bonded to a vapor-permeable polyurethane membrane, which acts as a barrier to liquid and microorganisms.

The silicone foam dressing include two types: silicone foam dressing with border and silicone foam dressing without border. The absorbent core of silicone foam dressing without border is polyurethane foam and that of silicone foam dressing with border consists of three components, a thin sheet of polyurethane foam ,a piece of nonwoven fabric and a layer of superabsorbent fibres. The core can lock exudate into the product, preventing the wound bed from drying out and the surrounding skin from becoming macerated.

Instructions

1. Clean the wound with saline
2. Dry skin around the wound, select the appropriate specification (crop the product if necessary). Cover the wound completely
3. Tear the face mask, stick the wound
4. Use the second dressing to fix (island)

Indications

The product is used for acute or chronic wounds and a variety of shades of wound.

1. Pressure ulcer
2. Venous leg ulcers
3. Diabetic foot ulcers
4. Burns 1st and 2nd degree
5. Skin graft donor site wound


Contraindications



1. The product is not suitable for dry wound
2. Infected wound must be used under control

Specification

TYPE	ITEM	SIZE	PCS/PKG	PKGS/CTN
Silicone Foam Dressing with Border	E7575	7.5*7.5cm	10	40
	E10100	10*10cm	10	40
	E12125	12.5*12.5cm	10	40
	E1515	15*15cm	10	40
	E1818	18*18cm	10	40
Silicone Foam Dressing without Border	F1010	10*10cm	10	40
	F1515	15*15cm	10	40

The product specifications can be available at clients' option.

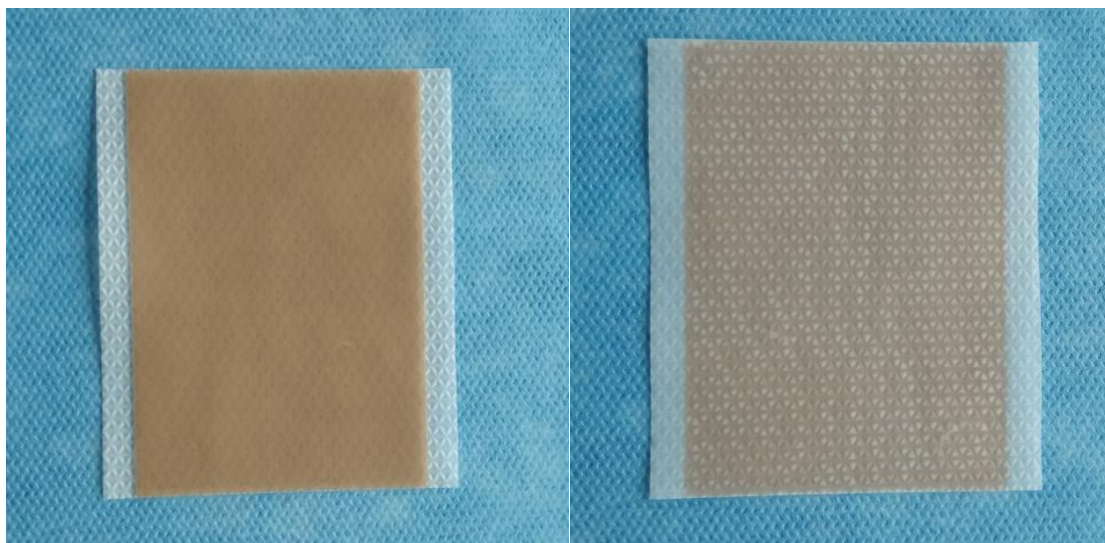

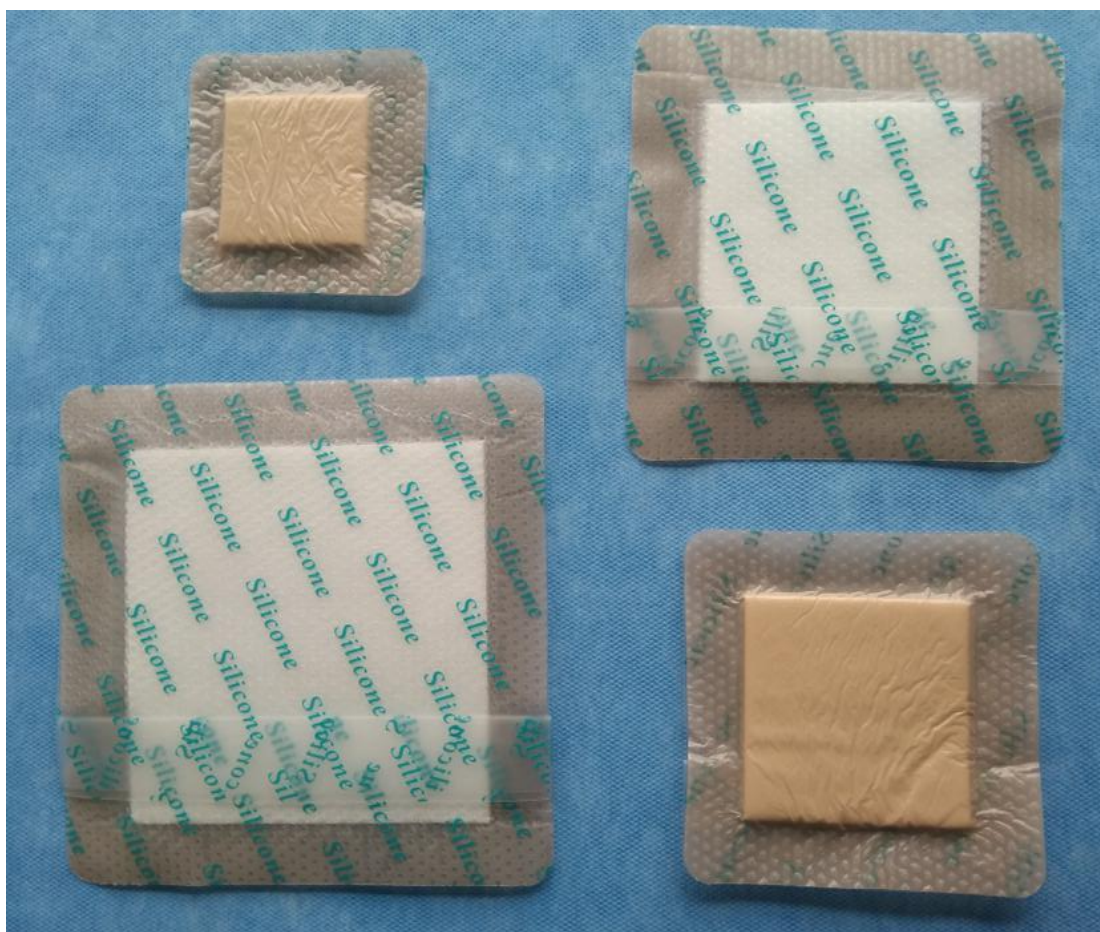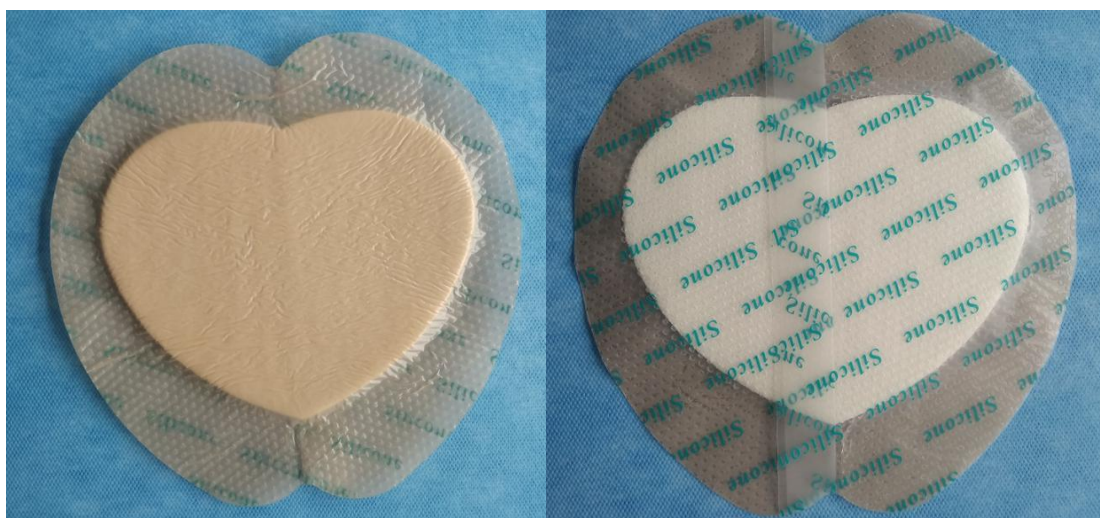

Barry Rong

Zibo Qichuang Medical Products Co., Ltd.

Mobile|Skype|Whatsapp|Wechat: +86 18560300608

Email: cn.barry@aliyun.com

Keywords: disposable nursing kit, catheter fixation device, medical dressing series, maternal and baby series, surgical tape, infusion paste, tourniquet